



Standard Features:

Construction

- Foot control with water on/off switch
- Grab handle for easy positioning
- Mount on Wall, Cabinet or Chair.
- Single pivoting instrument tray
- 3 full-size and 2 short modular electrical instrument panels for upgradability and serviceability
- Automatic handpiece holders
- Individual handpiece pressure and coolant adjustments
- Aircraft grade aluminum chassis construction
- Baked on Colorado Beige powder coat finish for durability
- Hospital grade plug and listed wiring and devices

Air Handpiece Connections

- Up to 3 non-optic 4-hole air handpieces

Syringes

- Air/water syringe
- Air-only Stropko™ syringe or air-only bonding syringe

Water System

- Dual closed water system mounts to wall, cabinet or chair
- DentaPure® water system disinfection – a single cartridge provides water treatment for one full year

Utilities

- Master air supply regulator and internal air filter
- 8' umbilical length, 1.25" outer diameter – with compressed air line and 1 main power cord

Custom Options:

Over-the-Patient Chair Adapter

Required for all radius mount chair units

- 2" diameter post
- A-Dec 500 adapter

Construction

- Plain disc foot control

Air Handpiece Connections

- Fiber optics (5-hole/6-pin) available – Halogen or LED options

Syringes

- Air/Water endodontic Stropko™ syringe
- 2" or 4" Stropko™ adapter

Endodontic Irrigation and Evacuation

- RootDown™ Irrigation – choose single, dual or triple irrigation
- RootDown™ Active Irrigation upgrade
- Micro evacuator

Assistant Package & Suction Connections

- Side or corner bar-mounted instrument holder bar
- High volume evacuator, straight, titanium
- Saliva ejector, titanium
- High volume evacuator, Ergo Angle, titanium
- Additional air/water syringe

Bracket Mounted Storage

- Vice bracket for Apex Locator – adds 6" to width
- Various instrument brackets (See below)


Computer Accessories

- USB port
- 15' USB cable

Utilities

- Floor junction box
- Extended foot control tubing

Illustrated Dimensions:

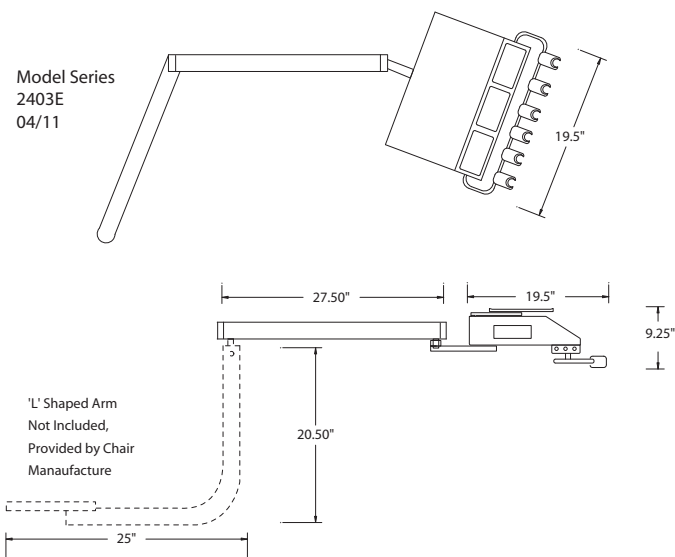

Dimensions with Optional Accessories:

Accessory Brackets Increase Width and/or Depth

- System B bracket – adds 2³/₈" to width
- Touch N Heat bracket – adds 1³/₈" to width
- Obtura gun bracket – adds 4¹/₄" to width
- Vice bracket for Apex Locator – adds 6" to width
- Irrigation (exterior mount) – adds 3.25" to width or depth