

When Precision Matters

DESIGNED FOR THE ORAL SURGEON
AND IMPLANT SPECIALIST



Proudly Made
in the USA!

ADVANCED
DENTAL
SYSTEMS®



ASI
Dental Specialties

844.880.3636 | ASIDental.com
8811 American Way, Suite 130 | Englewood, CO 80112

The Glider™

INTEGRATED DENTAL SURGICAL TABLE

Model 90-2148




An integrated instrument solution for efficient implant and surgical procedures. The Glider™ is easily positioned over standard and surgical style patient chairs and provides sufficient workspace to lay out sterile instruments while integrating surgical motors, handpieces, and suction instruments for ease of access by the dental team.

The Glider is mounted on stable, heavy-duty casters with a programmable electric height adjustment. This allows for a full range of operator capability for either sit-down or stand-up treatment with mobile positioning for desired placement.

STRUCTURAL FEATURES

- Smooth, easy-to-clean, non-porous solid surface worktop provides sufficient workspace for placing hand instruments and other devices
- Programmable electro-mechanically adjusted height with patient safety proximity sensor
- High-grade dual wheel casters for strength and ease of mobility
- Aluminum MetalCraft design for esthetics, durable lifetime quality and lighter weight
- Rear air and water quick connects and electrical outlets for table top instruments
- Rear-mounted holder for disposable suction canisters or solids trap filter for connection to central suction

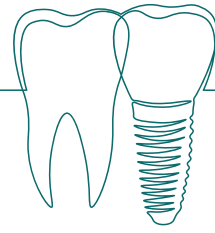

INSTRUMENT & ACCESSORY OPTIONS

- Integrated surgical motors with full range of handpiece attachments for surgical to implant
- Integrated piezo ultrasonic options
- Side-mounted dovetail bars on each side. Holds accessory items including saline hooks, sharps container, removable sterilizable handpiece holder, wipes bracket, and waste bins
- Side-mounted or rear mounted glove dispenser option
- Floating shelf to hold additional instruments or cassettes
- Heavy duty monitor mount or vitals monitor shelf

INSTRUMENT DELIVERY FEATURES

- Customizable instrument delivery with operation from only one foot control
- Standard pneumatic air handpiece connections and dual air/ water syringes for doctor and assistant
- Titanium suction handpieces with autoclavable tubing option
- Articulating handpiece holder bar for preferred placement
- Internal water bottle system for standard irrigation with surgical saline pump option



SCAN NOW FOR FULL
SPECIFICATIONS

844.880.3636 | ASIDental.com

OPTIONAL EQUIPMENT SOLUTIONS FOR THE GLIDER



VITALS MONITOR MOUNT

Swivel flat tray for holding vitals monitors at rear of Glider surgical table. Has basket to hold pressure cuff, pulse-ox, etc. Mounted on a vertical post riser. Platform measures 9" deep by 12" wide. [PN 90-1839]



HEAVY DUTY COMPUTER MONITOR MOUNT

Heavy duty monitor mount with articulating arm for varied positioning and height adjustment. [PN 90-1838]



DISPOSABLE SUCTION MOUNT

Adds a 1200cc disposable suction canister for use on the Glider. [PN 90-B040]



MOBILE UTILITY CART

Utility cart with internal vacuum and/or compressor to be used in office locations without central vacuum/compressor utility connections. Connects to Glider Surgical Table and plugs into standard electrical outlet for operation. [PN 90-9003]



SCAN NOW
FOR FULL
SPECIFICATIONS



844.880.3636 | ASIDental.com
8811 American Way, Suite 130 | Englewood, CO 80112